


## P<sup>3</sup> Functional Adult



## P<sup>3</sup> Functional Kiddythotics®

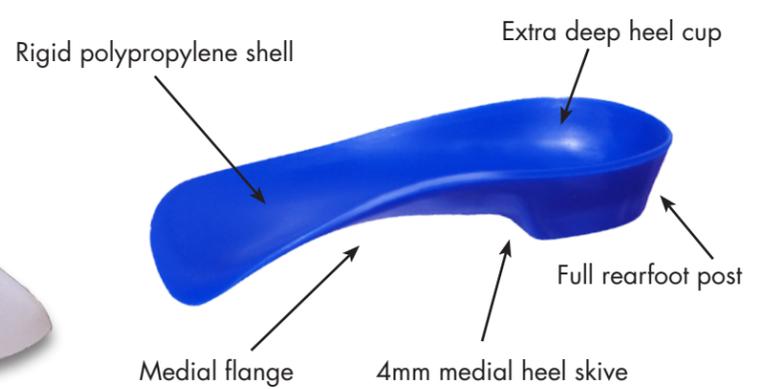

**PROLAB**  
ORTHOTICS

Evidence-Based  
Medicine

## Quick Guide P<sup>3</sup> Functional Prefabs



### Unposted (natural)

An excellent choice for shoes without a removable insole, patients with smaller feet, or casual shoes.



### Posted (black)

Select a posted prefab with a full rearfoot post for increased rearfoot and midfoot control.



### Kiddythotics®

(color-coded by size)  
Use Kiddythotics for pediatric flatfoot and calcaneal apophysitis.

### Polypropylene Shell

The most long-lasting, durable material that resists compression. Ensures consistent control without flattening.

Semirigid

Semirigid

Rigid

### Medial Heel Skive

Supports the medial column to attenuate pronation.

2 mm

2 mm

4 mm

### Rearfoot Post

Increases rearfoot and midfoot stability.

Spot Grind

Full

Full

### Deep Heel Cup

Improves calcaneal stabilization.

14 mm

14 mm

Extra Deep

### Forefoot Valgus Correction

Stabilizes the midtarsal joint to reduce forefoot motion addressing the valgus forefoot deformity found in most pathologic feet.

Included

Included

Not Available

### Medial Flange

Holds foot in healthy position for arch development. May prevent long-term problems due to pronation.

Not Available

Not Available

Included

## Options

### Topcover

Full-length topcovers provide cushioning and comfort, and allow you to add forefoot corrections. P<sup>3</sup> Prefabs with topcovers are ideal for replacing removable insoles in athletic shoes.



### Unposted TC

Full-length tricolor EVA topcover (TC) on an unposted P<sup>3</sup> Prefab provides extra cushioning.

### Posted TC

Full-length tricolor EVA topcover (TC) on a posted P<sup>3</sup> prefab provides extra cushioning.

### Sport

Includes a full-length, high-rebound EVA and microsuede topcover on a P<sup>3</sup> posted prefab. Ideal for active patients and athletes.

### Low-Impact

Includes a full-length plastazote and Poron™ topcover on a P<sup>3</sup> posted prefab for patients with sensitive feet needing extra cushioning.

### Kiddythotics® Sizer Set

Accurate sizing is easy using the Kiddythotics sizer set. Includes one prefab of each size.

