


P³ Functional Adult



P³ Functional Kiddythotics®

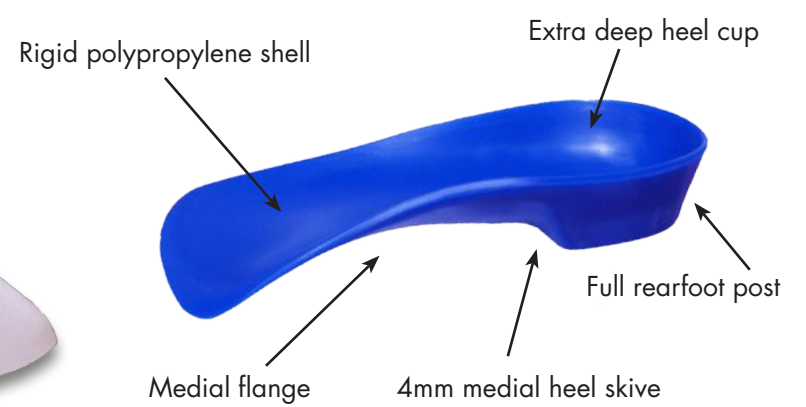

PROLAB
ORTHOTICS

Evidence-Based
Medicine

Quick Guide P³ Functional Prefabs



Unposted (natural)

An excellent choice for shoes without a removable insole, patients with smaller feet, or casual shoes.



Posted (black)

Select a posted prefab with a full rearfoot post for increased rearfoot and midfoot control.



Kiddythotics® (color-coded by size)

Use Kiddythotics for pediatric flatfoot and calcaneal apophysitis.

Polypropylene Shell

The most long-lasting, durable material that resists compression. Ensures consistent control without flattening.

Semirigid

Semirigid

Rigid

Medial Heel Skive

Supports the medial column to attenuate pronation.

2 mm

2 mm

4 mm

Rearfoot Post

Increases rearfoot and midfoot stability.

Spot Grind

Full

Full

Deep Heel Cup

Improves calcaneal stabilization.

14 mm

14 mm

Extra Deep

Forefoot Valgus Correction

Stabilizes the midtarsal joint to reduce forefoot motion addressing the valgus forefoot deformity found in most pathologic feet.

Included

Included

Not Available

Medial Flange

Holds foot in healthy position for arch development. May prevent long-term problems due to pronation.

Not Available


Not Available

Included

Options

Topcover

Full-length topcovers provide cushioning and comfort, and allow you to add forefoot corrections. P³ Prefabs with topcovers are ideal for replacing removable insoles in athletic shoes.

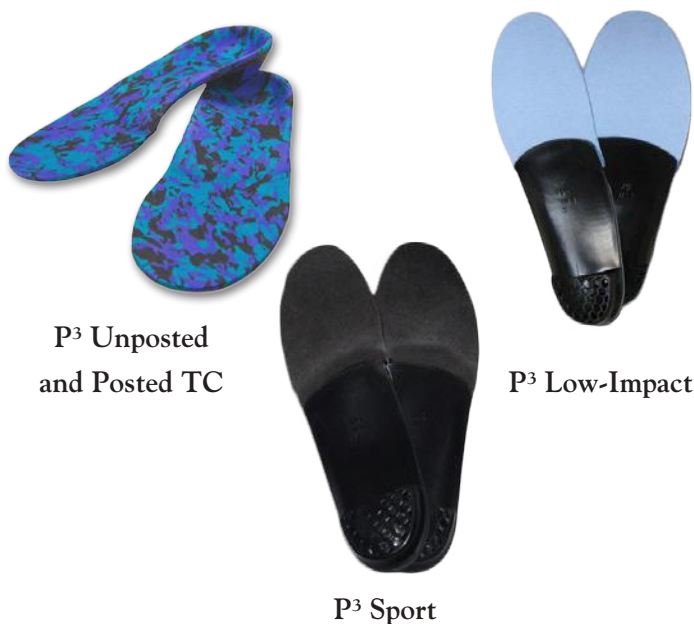

Unposted TC

Full-length tricolor EVA topcover (TC) on an unposted P³ Prefab provides extra cushioning.

Posted TC

Full-length tricolor EVA topcover (TC) on a posted P³ prefab provides extra cushioning.

Sport

Includes a full-length, high-rebound EVA and microsuede topcover on a P³ posted prefab. Ideal for active patients and athletes.

Low-Impact

Includes a full-length plastazote and Poron™ topcover on a P³ posted prefab for patients with sensitive feet needing extra cushioning.

Kiddythotics® Sizer Set

Accurate sizing is easy using the Kiddythotics sizer set. Includes one prefab of each size.

